

NOVO Go⁺

User Manuals

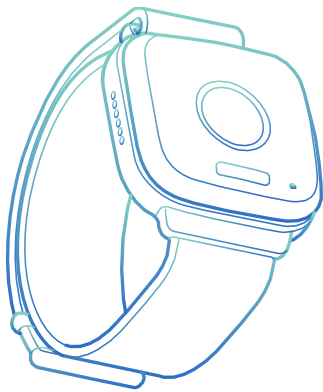
NOVO Go

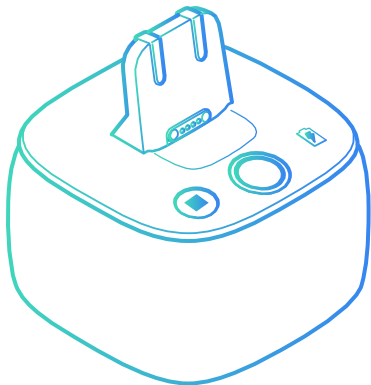NOVO Go DOCK

[English](#)

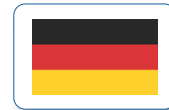[Deutsch](#)

[Español](#)

[Svenska](#)


[Français](#)


English

User Manual

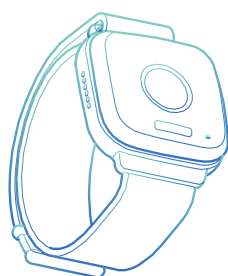[NOVO Go](#)

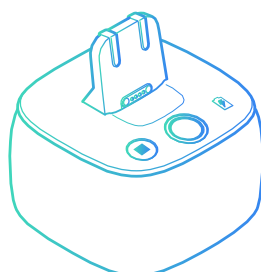
[NOVO Go DOCK](#)

NOVO GO^{III}

Portable Carephone

User Manual

NE41 21002-02 v2.0

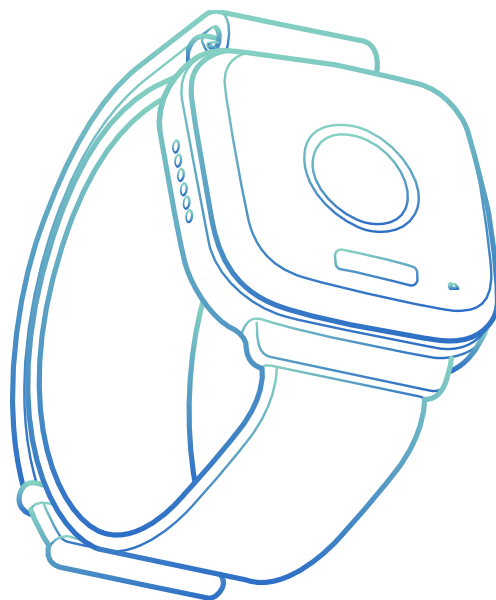Declaration of Conformity: Hereby Legrand AB declares that this device is in compliance with Directive 2014/53/EU. The full text of the EU Declaration of Conformity is available at:
www.neat-group.com/downloads/declarations-of-conformity

Table of Contents

1 About NOVO Go	5
1.1 Contents in the Box	6
2 Your NOVO Go	7
3 Get Started with NOVO Go	8
3.1 Get started with NOVO Go + NOVO Go DOCK	8
4 Accessory	9
4.1 Attach Accessory	9
5 Send Alarm with NOVO Go	10
5.1 Send Alarm with Pre-Call	10
5.2 Cancel Alarm	11
5.3 Send Alarm without Pre-Call	11
6 Incoming Voice Call	12
6.1 Answer Voice Call Manually	12
6.2 Reject Incoming Voice Call	12
7 Charge NOVO Go	13
7.1 Charge with USB Charger	13
8 Operational Functions	15
8.1 Standby Mode	15
8.2 Status Check	15
8.3 Mobile Network Connection Status	16
8.4 Device Temperature	16
8.5 Set NOVO Go to On	16
8.6 Set NOVO Go to Off	16
9 More Functions	18
9.1 Passive Alarm (Optional)	18
9.2 I'm OK (Optional)	18
9.3 Reminder (Optional)	19
9.4 Find My NOVO Go	19
10 Safety and Handling Information	20
10.1 Cleaning	20

10.2 Discard	21
11 Technical Data	22
Appendix 1 NOVO Go + Legrand Care Stationary Carephone	23
Appendix 2 Connect with Alarm Devices	24
1 Send Alarm with Alarm Device	24
2 Answer Voice Call with Alarm Device	24

1

About NOVO Go


NOVO Go is a portable carephone with integrated SIM card and GPS that gives information about the position of the device. Send an alarm with the easy-to-press button and have a two-way voice communication with the alarm receiver. You can wear NOVO Go with a wristband or as a necklace.

Use NOVO Go as follows:

- stand-alone carephone.
- together with NOVO Go DOCK.
- together with a Legrand Care stationary carephone.


For information about NOVO Go DOCK and the Legrand Care stationary carephones, refer to their related user manuals.

NOVO Go has radio extension to other devices and is fully integrated with Pulse CMP and the IoT Telecare ecosystem supplied by Legrand Care.

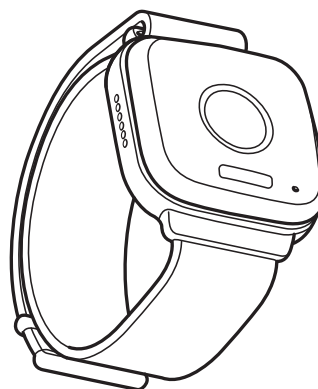This manual describes the basic usage of NOVO Go. For more information on how to do adjustments and configurations, speak with your Legrand Care distributor.

1.1 Contents in the Box

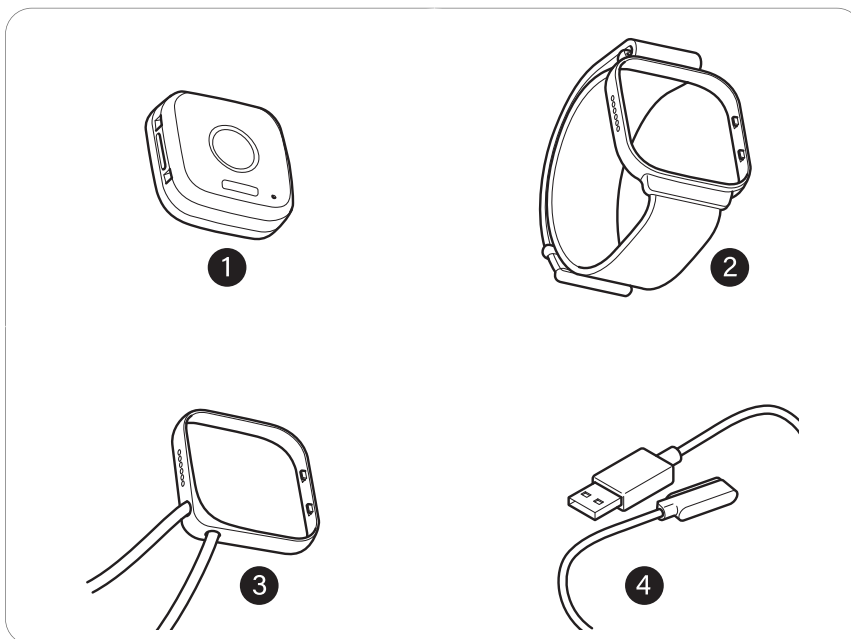
- 1 NOVO Go Device
- 2 Wristband Accessory
- 3 Necklace Accessory
- 4 USB Charger (included AC Power Adapter for EU or UK)

NOVO Go Quick Reference Guide

2

Your NOVO Go

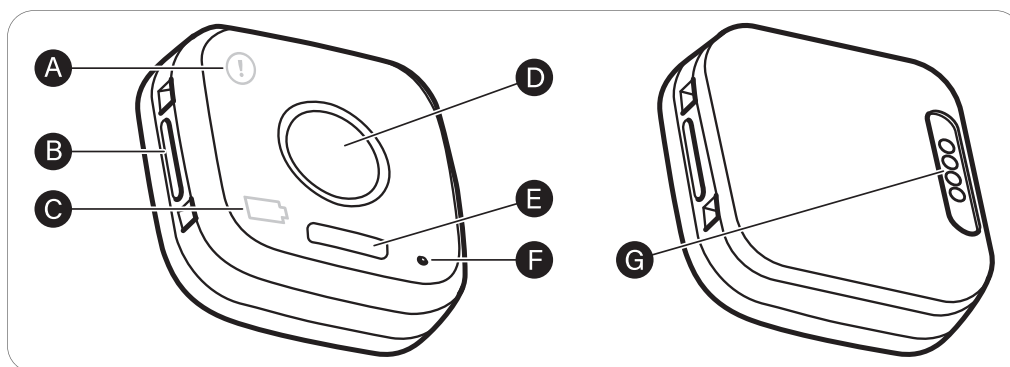- A** Network Icon
- B** Loudspeaker
- C** Battery Icon
- D** Alarm Button
- E** Extra Button
- F** Microphone
- G** Charge Connector

RELATED CHAPTER(S)

[5 Send Alarm with NOVO Go](#)

[8 Operational Functions](#)

3

Get Started with NOVO Go

- 1 Remove NOVO Go from the box.
- 2 Examine the contents for defects or missing items. If something is missing or if any item has defects, speak with your Legrand Care distributor.
- 3 Attach an accessory to NOVO Go, refer to [4 Accessory](#).
- 4 Charge NOVO Go, refer to [7.1 Charge with USB Charger](#).
- 5 NOVO Go starts automatically. If it does not start, refer to [8.5 Set NOVO Go to On](#).
- 6 The software updates automatically. The Alarm button flashes yellow and the Network icon flashes green for up to 10 minutes.
- 7 After the Alarm button goes off, continue with step 8.
- 8 Do a status check of NOVO Go, refer to [8.2 Status Check](#).
- 9 Do an operational test of NOVO Go, refer to [5 Send Alarm with NOVO Go](#).
- 10 If the test result is unsatisfactory, speak with your Legrand Care distributor.
- 11 If the test result is satisfactory, you are good to go!

3.1 Get started with NOVO Go + NOVO Go DOCK

Refer to NOVO Go DOCK User Manual.


4

Accessory

The same procedure is applicable to the two accessories.

4.1 Attach Accessory

- 1 Put NOVO Go on something stable, with the alarm button side up.
- 2 Position the accessory above NOVO Go. Refer to the illustration below.
- 3 Make sure that the holes on the side **1** point in the direction of the loudspeaker **B**.
Refer to [2 Your NOVO Go](#)
- 4 Push until you hear click sounds on the two sides.

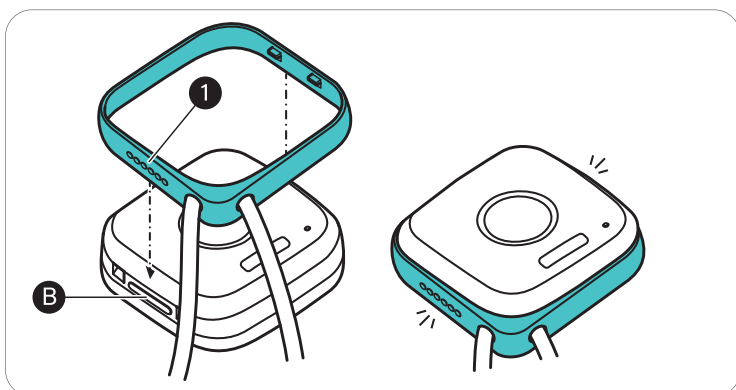RELATED CHAPTER(S)

[3 Get Started with NOVO Go](#)


5

Send Alarm with NOVO Go


In the basic configuration NOVO Go sends an alarm with a pre-call. To change the pre-call configuration, speak to your Legrand Care distributor.

5.1 Send Alarm with Pre-Call


In this configuration an incorrect alarm can be cancelled, refer to [5.2 Cancel Alarm](#).

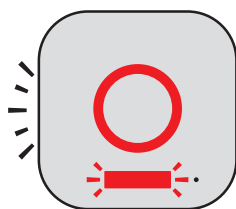- 1 Press the Alarm button **D**.
 - The Alarm button goes red.
 - The Extra button **E** flashes red.
 - A pre-call alarm operates.
 - It is possible to cancel the alarm.

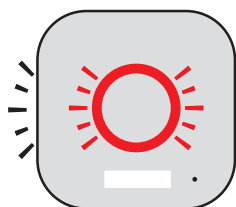- 2 After the pre-call alarm NOVO Go sends the alarm.
 - The Alarm button flashes red.
 - NOVO Go makes a sound.

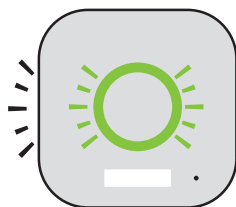- 3 When the alarm is received, the Alarm button changes to green.
- 4 Speak with the person who receives your alarm. The call ends when the person who receives the alarm hangs up.

RELATED CHAPTER(S)

[2 Your NOVO Go](#)

[3 Get Started with NOVO Go](#)

5.2 Cancel Alarm


During the pre-call alarm, press the Extra button **E**. The sound stops and the Extra button and the Alarm button go off.

RELATED CHAPTER(S)

[2 Your NOVO Go](#)

[5 Send Alarm with NOVO Go](#)

5.3 Send Alarm without Pre-Call

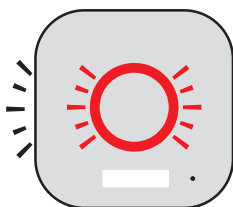- 1 Press the Alarm button **D**.
- 2 NOVO Go makes a short sound and sends the alarm.
- 3 During the transmission of the alarm:
 - The Alarm button flashes red.
 - NOVO Go makes a sound.

- 4 When the alarm is received, the Alarm button changes to green.
- 5 Speak with the person who receives your alarm. The call ends when the person who receives the alarm hangs up.

RELATED CHAPTER(S)

[2 Your NOVO Go](#)

[3 Get Started with NOVO Go](#)


6

Incoming Voice Call


NOVO Go can receive voice calls from approved phone numbers.

6.1 Answer Voice Call Manually

In this configuration an incoming voice call can be rejected, refer to [6.2 Reject Incoming Voice Call](#).

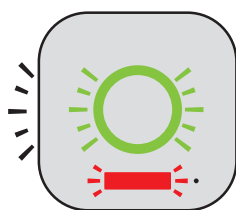1 NOVO Go gives a call signal. The Alarm button flashes green and the Extra button flashes red.

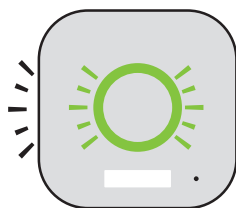
2 Press the Alarm button. The Alarm button stays green and the Extra button goes off.

3 Speak with the person who calls. The voice call ends when the person who called hangs up, or if the extra button on NOVO Go is pressed.

6.2 Reject Incoming Voice Call

Reject an incoming voice call as follows:

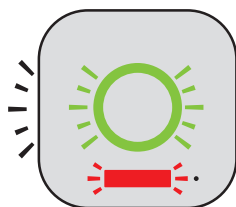1 NOVO Go gives a call signal. The Alarm button flashes green and the Extra button flashes red.

2 Press the Extra button.

- The call signal stops.
- The Alarm button and the Extra button go off.

7


Charge NOVO Go

If NOVO Go is set to off (for example if the battery has not been charged), it automatically sets to on when you connect to a charger. If the battery icon slowly flashes red, it is necessary to charge the battery.

7.1 Charge with USB Charger

- 1 Connect the USB cable to the AC power adapter and attach the magnetic charger to the back of NOVO Go.

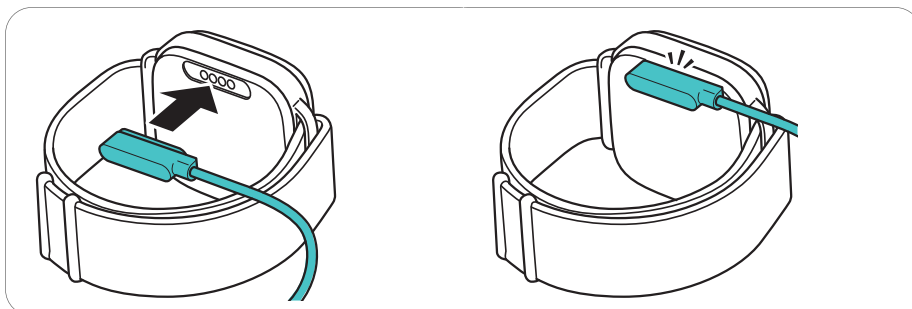
- 2 Connect the AC power adapter to the wall socket. NOVO Go gives a short sound. The battery icon flashes green until the battery is fully charged.

- 3 If necessary, do a check of the status during the charging, refer to [8.2 Status Check](#).
- 4 When NOVO Go is fully charged, the battery icon stays on for 1 minute and then it goes off.

> 1 minute >

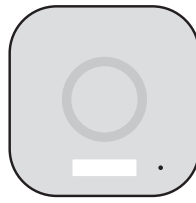RELATED CHAPTER(S)

[3 Get Started with NOVO Go](#)

8

Operational Functions

8.1 Standby Mode

When NOVO Go is in standby mode, all lights are set to Off.

After a procedure, NOVO Go goes back to standby mode after 2 to 5 seconds.

8.2 Status Check

Do a status check of NOVO Go as follows:

- 1 Press the Extra button **E**.
- 2 The Alarm button shows one condition as follows:

Red:
Battery charged
21-40%

Yellow:
Battery charged
41-60%

Green:
Battery charged
61-100%

RELATED CHAPTER(S)

[2 Your NOVO Go](#)

[3 Get Started with NOVO Go](#)

[7.1 Charge with USB Charger](#)

8.3 Mobile Network Connection Status


If the Network icon flashes, NOVO Go cannot find connection to the mobile network and no alarm can be sent.

The Network icon goes off when NOVO Go finds connection.

8.4 Device Temperature

If NOVO Go gets too hot, all LED lights start to flash and most functions go off. Correct functionality is restored when NOVO Go gets a lower temperature.

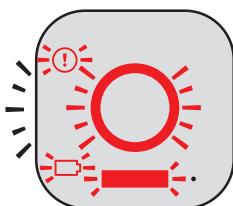8.5 Set NOVO Go to On

- 1 Press and hold the Alarm button **D** and the Extra button **E** at the same time for 5 seconds.
- 2 The Alarm button comes on with a green light for a period of 2 seconds.
- 3 The Network icon slowly flashes for a period of 5 seconds, refer to [8.3 Mobile Network Connection Status](#)

8.6 Set NOVO Go to Off

- 1 Press and hold the Alarm button **D** and the Extra button **E** at the same time. NOVO Go makes a sound.
- 2 After 5 seconds the Alarm button comes on with a red light and then slowly goes off.

**NOTE**

If you set NOVO Go to off, it does not receive radiosignals from connected devices, for example smoke detectors and NOVO Go DOCK.


9

More Functions

9.1 Passive Alarm (Optional)

If NOVO Go has not sensed any movement during a period of time (usually for 24 hours), it operates as follows:

- 1 The Extra button flashes yellow for 60 seconds and NOVO Go tells you to press the Extra button.

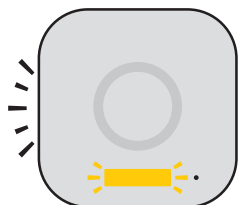
- 2 Press the Extra button. The alarm cancels.
- 3 If you do not press the Extra button, NOVO Go sends an alarm.

9.2 I'm OK (Optional)

In a configuration with the "I'm OK" function you must press the Extra button once a day during a specified time period, in order to stop an alarm to be sent from NOVO Go.

- 1 During all of the specified "I'm OK" time period: The Extra button flashes yellow once per second.
- 2 During the last 30 minutes of the "I'm OK" period: The Extra button flashes and NOVO Go gives a short sound signal every 5 minutes.
- 3 During the last 5 minutes of the "I'm OK" period: The Extra button flashes, NOVO Go gives a short sound signal once a minute and says "Please press the Extra button".

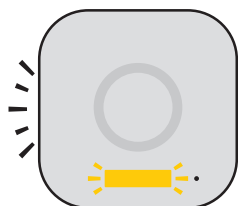- 4 Press the Extra button. The "I'm OK" function is reset.
- 5 If you do not press the Extra button, NOVO Go sends an alarm.

9.3 Reminder (Optional)

In this configuration NOVO Go tells you to do a task (for example to take medication). NOVO Go operates as follows:

- 1 A voice tells you a task. The Extra button flashes green.

- 2 Do the requested task.
- 3 Press the Extra button. The voice stops and the Extra button goes off.

9.4 Find My NOVO Go

This configuration is only available together with NOVO Go DOCK, refer to the NOVO Go DOCK User Manual.

10

Safety and Handling Information

- Read the manual before use.
- Do not use the device if there is a risk of death, call Emergency number for applicable medical care.
- Do an operational test before use and after adjustments.
- Only use the device in temperatures between -20 to +50 °C. Avoid exposure to direct sunlight and do not use near fire.
- Only connect with other devices, chargers and power adapters that are approved by Legrand Care.
- Set the device to off if it will not be used for a longer period of time, refer to [8.6 Set NOVO Go to Off](#).
- Be careful with the device. Do not hit, shake, twist or let it fall.
- Make sure that the USB cable and the AC power cable cannot be stepped on, tripped over or otherwise subjected to damage or stress.
- Do not hold the device close to your ears or to your hearing aid. The sound can be high and cause damage to you ears.
- Do not open the device. For repairs, speak with your Legrand Care distributor.
- Do not put NOVO Go and/or the USB charger near items that are sensitive to magnets, for example credit cards, bank cards and memory devices.
- This device contains a connection to the mobile network and should be treated in the same way as a mobile phone, for example in an airplane.
- To give the accurate location, the GPS signal needs to be available. The precision of the location is limited indoors.
- If you have a pacemaker or other implanted medical device, speak with your physician before you use NOVO Go.

10.1 Cleaning

- Use a weak soap solution or a surface disinfection and a moist cloth to clean the device.
- Do not use strong chemicals, grease and other harsh material.
- After cleaning, do an operational test, refer to [5 Send Alarm with NOVO Go](#)

10.2 Discard

Do not put the used device and/or batteries in the household waste. If you put the device and/or batteries in the household waste, they can cause damage to the environment. To discard and recycle the device and/or batteries, refer to local regulations or your Legrand Care distributor.

11

Technical Data

Measures (W × H × D)	40 × 40 × 13 mm
Weight	33 g
Operating Temperature	Between -20 and +50 °C
Battery	Lithium Polymer, 370 mAh
Battery Charging Time	Up to 4.5 hours
Communication	2G/3G/4G, VoLTE, roaming SIM
Positioning	Assisted/standalone GPS, Galileo, GLONASS, BeiDou, QZSS
GPS Precision	~5 meters open sky
Frequency _{RF} (EU/UK)	869.200–869.250 MHz (Social Alarms) Category 1
Transmit Power _{max} (EU/UK)	869.2 MHz – 10 mW (EIRP)
IP Class	IP54
Power Supply	Magnetic USB connector, 5 Vdc, 1A
Supply Voltage	5 V
Battery performance	~24 hours in normal use

Appendix 1. NOVO Go + Legrand Care Stationary Carephone

It is possible to use NOVO Go as a pendant to a Legrand Care stationary carephone when the user is at home. Contact your Legrand Care distributor.

Appendix 2. Connect with Alarm Devices

NOVO Go can be connected with other alarm devices, for example SMILE, TOUCH, DOOR and PIR II. Contact your Legrand Care distributor.

1 Send Alarm with Alarm Device

On the connected alarm device, press the alarm button. NOVO Go sends the alarm, with position details, to the alarm central, refer to [5 Send Alarm with NOVO Go](#).

2 Answer Voice Call with Alarm Device

To answer an incoming voice call, press the alarm button on the connected alarm device.

Refer to [6 Incoming Voice Call](#).

All rights reserved. Circle font is used with attribution to <https://www.photoshopsupply.com>

Document ID: NE41 21002-02 v2.0

Revision date: 2023-12-08

Legrand AB

Varuvägen 2

246 42 Löddeköpinge

+46 46 707065

www.legrandcare.com

infosweden.legrandcare@legrand.com

NOVO Go⁺

NOVO Go DOCK

User Manual

NE41 21003-02 v2.0

↑
↑

[Top](#)

↑

[Top - English](#)

 **legrand**[®] | care

Declaration of Conformity: Hereby Legrand AB declares that this device is in compliance with Directive 2014/53/EU. The full text of the EU Declaration of Conformity is available at:
www.neat-group.com/downloads/declarations-of-conformity

Table of Contents

1	About NOVO Go DOCK	4
1.1	Contents in the Box	4
2	Your NOVO Go DOCK	5
3	Get started with NOVO Go DOCK	6
3.1	Set NOVO Go DOCK to ON	7
3.2	Set NOVO Go DOCK to Off	7
3.3	Pair NOVO Go DOCK and NOVO Go	8
4	How to use NOVO Go DOCK	9
4.1	Charge NOVO Go with NOVO Go DOCK	9
4.2	NOVO Go DOCK as Loudspeaker	9
4.3	Send Alarm with NOVO Go DOCK	10
4.4	Find My NOVO Go	10
4.5	Backup Battery	11
5	Safety and Handling Information	12
5.1	Cleaning	12
5.2	Discard	13
6	Technical Data	14

1

About NOVO Go DOCK


NOVO Go DOCK is a device that gives NOVO Go more functions as follows:

- Easy-to-use charging.
- Loudspeaker with a loud and clear sound when at home.
- Backup battery that charges NOVO Go (for example if a power failure occurs).
- Auto home/away function to keep the battery in NOVO Go charged for a longer period of time.
- Find My NOVO Go function.

NOVO Go DOCK does not have its own connectivity and must be used with a NOVO Go to send an alarm.

This user manual gives information about the basic configuration of NOVO Go DOCK and how to use it together with NOVO Go. For more information on how to set up and do configurations, speak with your Legrand Care distributor.


1.1 Contents in the Box

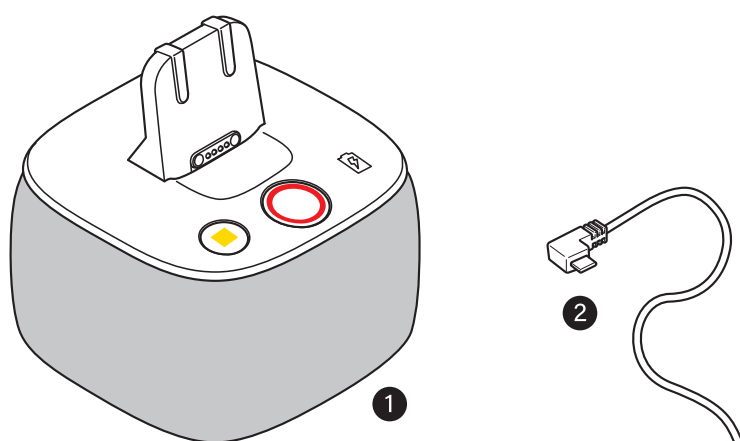- 1 NOVO Go DOCK
- 2 Micro-USB with AC Power Adapter for EU or UK

2

Your NOVO Go DOCK

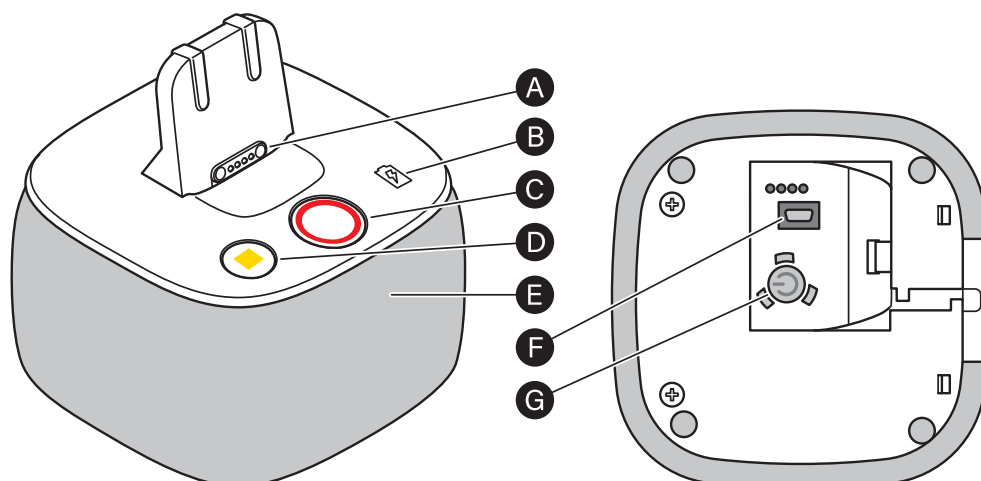- A** Charge Connector
- B** Power Indicator
- C** Alarm Button
- D** Extra Button
- E** Loudspeaker
- F** Micro-USB connection (on the bottom, under the lid)
- G** On/Off Button (on the bottom, under the lid)

RELATED CHAPTER(S)

[3.1 Set NOVO Go DOCK to ON](#)

[3.2 Set NOVO Go DOCK to Off](#)

[3.3 Pair NOVO Go DOCK and NOVO Go](#)


[4.3 Send Alarm with NOVO Go DOCK](#)

[4.4 Find My NOVO Go](#)


3

Get started with NOVO Go DOCK

- 1 Remove NOVO Go DOCK from the box.
- 2 Examine the contents for missing or defect items. If something is missing or defect, speak with your Legrand Care distributor.
- 3 Keep NOVO Go DOCK on something stable near a wall socket.
- 4 Open the lid on the bottom of NOVO Go DOCK.

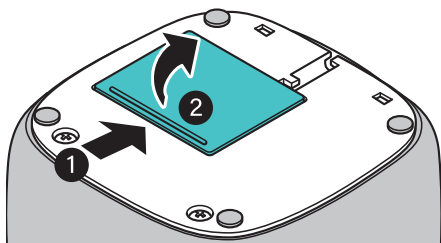- 5 Remove the plastic protection from the Micro-USB connector.
- 6 Connect the Micro-USB connector and carefully attach the cable to the hooks, refer to the image below.

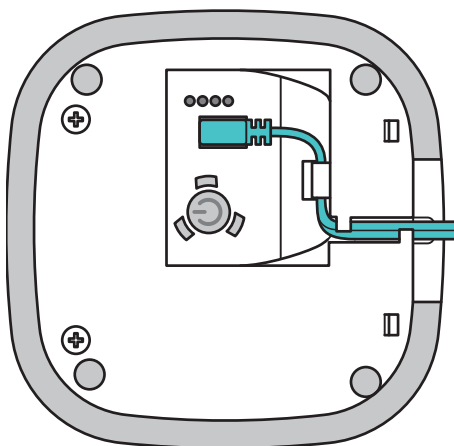- 7 Connect the AC adapter to the wall socket. NOVO Go DOCK sets to on automatically.
- 8 Pair NOVO Go DOCK with NOVO Go, refer to [3.3 Pair NOVO Go DOCK and NOVO Go](#).
- 9 Put NOVO Go in NOVO Go DOCK.

- 10 Do a test to make sure that NOVO Go DOCK operates correctly, refer to [4.3 Send Alarm with NOVO Go DOCK](#).
- 11 Remove NOVO Go from NOVO Go DOCK and do step 10 again.
- 12 If NOVO Go does not send alarms, speak with your Legrand Care distributor.
- 13 If NOVO Go sends alarms, you are good to go!

3.1 Set NOVO Go DOCK to ON

- 1 Open the lid on the bottom of NOVO Go DOCK.
- 2 Press and hold the On/Off button  for more than 3 seconds. The buttons on the top of NOVO Go DOCK flash and then go off.
- 3 After approximately 3 seconds the Alarm button comes on and stays red.

NOTE

This applies to standard configuration.

- 4 Close the lid.

RELATED CHAPTER(S)

[2 Your NOVO Go DOCK](#)

[3 Get started with NOVO Go DOCK](#)

3.2 Set NOVO Go DOCK to Off

- 1 Open the lid on the bottom of NOVO Go DOCK.
- 2 Press and hold the On/Off button .
- 3 After 3 seconds NOVO Go DOCK makes a sound and the Alarm button goes off.
- 4 Close the lid.


RELATED CHAPTER(S)

[2 Your NOVO Go DOCK](#)

3.3 Pair NOVO Go DOCK and NOVO Go

Pair NOVO Go DOCK and NOVO Go manually as follows:

- 1 Make sure that the holes on the side **1** point up.
- 2 Put NOVO Go in NOVO Go DOCK.

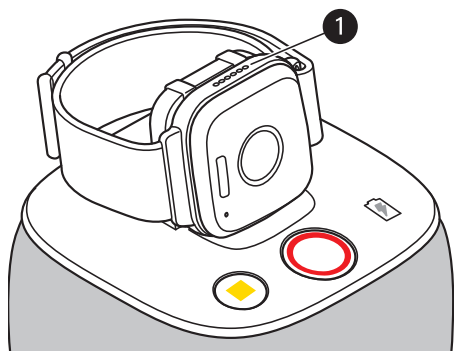- 3 NOVO Go makes a sound and the Battery icon flashes green.

- 4 On NOVO Go DOCK, press and hold the Alarm button **C** and the Extra button **D** at the same time for more than 3 seconds.
- 5 NOVO Go DOCK makes a sound, and is paired with NOVO Go.

RELATED CHAPTER(S)

[2 Your NOVO Go DOCK](#)

[3 Get started with NOVO Go DOCK](#)


[4.2 NOVO Go DOCK as Loudspeaker](#)

4


How to use NOVO Go DOCK

4.1 Charge NOVO Go with NOVO Go DOCK

- 1 Make sure that the holes on the side ❶ point up.
- 2 Put NOVO Go in NOVO Go DOCK.

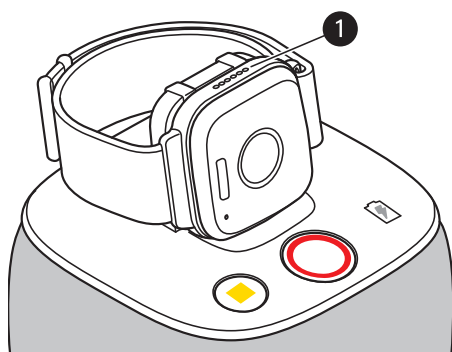- 3 NOVO Go makes a short sound and the Battery icon flashes green.

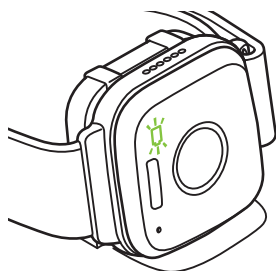- 4 When NOVO Go is fully charged the Battery icon will stay on for 1 minute and then go off.

4.2 NOVO Go DOCK as Loudspeaker

Put NOVO Go in NOVO Go DOCK to make NOVO Go DOCK operate as a loudspeaker.

4.3 Send Alarm with NOVO Go DOCK

NOTE

Before you use this function for the first time, make sure that NOVO Go DOCK and NOVO Go are paired, refer to [3.3 Pair NOVO Go DOCK and NOVO Go](#).

Use NOVO Go DOCK to send an alarm as follows:

- 1 Make sure that NOVO Go DOCK and NOVO Go are in the same location. It is not necessary to put NOVO Go in NOVO Go DOCK to use this function.
- 2 On NOVO Go DOCK, press the Alarm button **C**. NOVO Go DOCK tells NOVO Go to send the alarm.
- 3 NOVO Go sends the alarm. Refer to the configuration of your NOVO Go.

RELATED CHAPTER(S)

[2 Your NOVO Go DOCK](#)


[3 Get started with NOVO Go DOCK](#)

[5 Safety and Handling Information](#)

4.4 Find My NOVO Go

If your NOVO Go is missing, find it as follows:

- 1 Press and hold the Extra button **D** for 3 seconds. The Extra button flashes yellow.
- 2 The missing NOVO Go makes a clear sound at regular intervals. The Alarm button and the Extra button flash green during a period of 60 seconds.

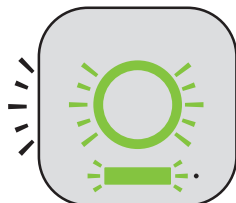- 3 If the missing NOVO Go is found, press the Extra button on your NOVO Go. The sound stops and the Extra button and Alarm button go off.
- 4 The Extra button on NOVO Go DOCK goes off automatically after 60 seconds.

RELATED CHAPTER(S)

[2 Your NOVO Go DOCK](#)

4.5 Backup Battery

If necessary (for example if a power failure occurs or if the AC adapter is disconnected), NOVO Go DOCK operates as follows:

- Continues to charge NOVO Go for approximately 72 hours.
- Sets the Alarm button lighting on NOVO Go DOCK to off to save power.

The Power indicator shows the status of the backup battery as follows:

Power available. The power indicator is set to off. The backup battery operates correctly.

No power. The power indicator slowly flashes green. The backup battery operates correctly. Make sure that the AC adapter is connected to the wall socket.

No power. The power indicator slowly flashes red. The power in the backup battery is low. Make sure that the AC adapter is connected to the wall socket.

Power available. The power indicator comes on with a red light. There is a malfunction of the backup battery. Keep the AC adapter connected to the wall socket and speak with your Legrand Care distributor.

5

Safety and Handling Information

- Read the manual before use.
- Do not use the device if there is a risk of death, call Emergency number for applicable medical care.
- Do an operational test before use and after adjustments.
- Only use the device in temperatures between 0 and +40 °C. Avoid exposure to direct sunlight and do not use near fire.
- Only connect with other devices, chargers and power adapters that are approved by Legrand Care.
- Keep away from dust, moist and dirt.
- Be careful with the device. Do not hit, shake, twist or let it fall.
- Make sure that the USB cable and the AC power cable cannot be stepped on, tripped over or otherwise subjected to damage or stress.
- Do not hold the device close to your ears or to your hearing aid. The sound can be high and cause damage to you ears.
- For repairs, speak with your Legrand Care distributor.

5.1 Cleaning

- Use a weak soap solution or a surface disinfection and a moist cloth to clean the device.
- Do not use strong chemicals, grease and other harsh material.

5.2 Discard

Do not put the used device and/or batteries in the household waste. If you put the device and/or batteries in the household waste, they can cause damage to the environment. To discard and recycle the device and/or batteries, refer to local regulations or your Legrand Care distributor.

6

Technical Data

Measures (W × H × D)	100 × 100 × 105 mm
Weight	440 g
Operating Temperature	Between 0 and +40 °C
Battery	NiMH, 2000mAh
Frequency _{RF} (EU/UK)	869.2 MHz (Social Alarms) Category 1
Transmit Power _{max} (EU/UK)	10 mW (EIRP)
Power Supply	5 Vdc, 1.0-1.5 A
Supply Voltage	5 V

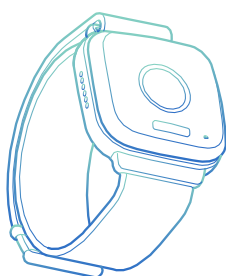
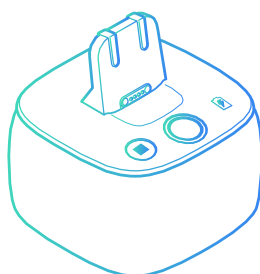
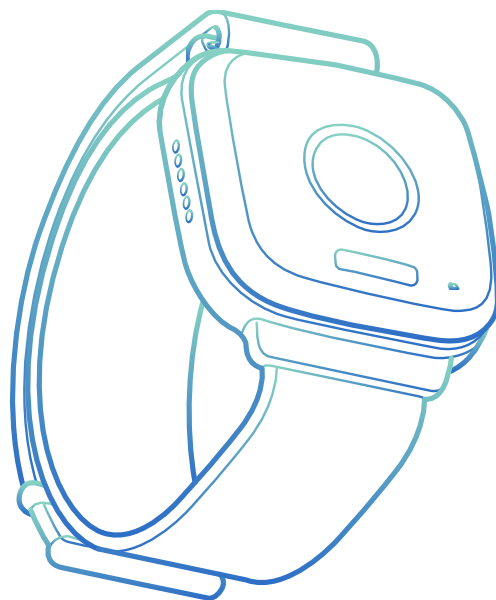
All rights reserved. Circle font is used with attribution to <https://www.photoshopsupply.com>

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